


# HD balance

## Next generation tilt-in-space wheelchairs!

HD Balance is a brand new Swedish tilt-in-space wheelchair with a design in which both users' and carers' needs and requirements have been considered throughout the developing process. Exclusively materials and components of the highest quality standards are used in the manufacturing of the wheelchair. Precision and durability are key words.

Seat widths: 38, 42, 46, 50  
Wheel sizes: 16" & 24"  
Seat tilt: ca 0° - 20° (standard)  
ca -5° - 30° (optional)  
Backrest angle: ca 90° - 120°  
Max user weight: 135 kg

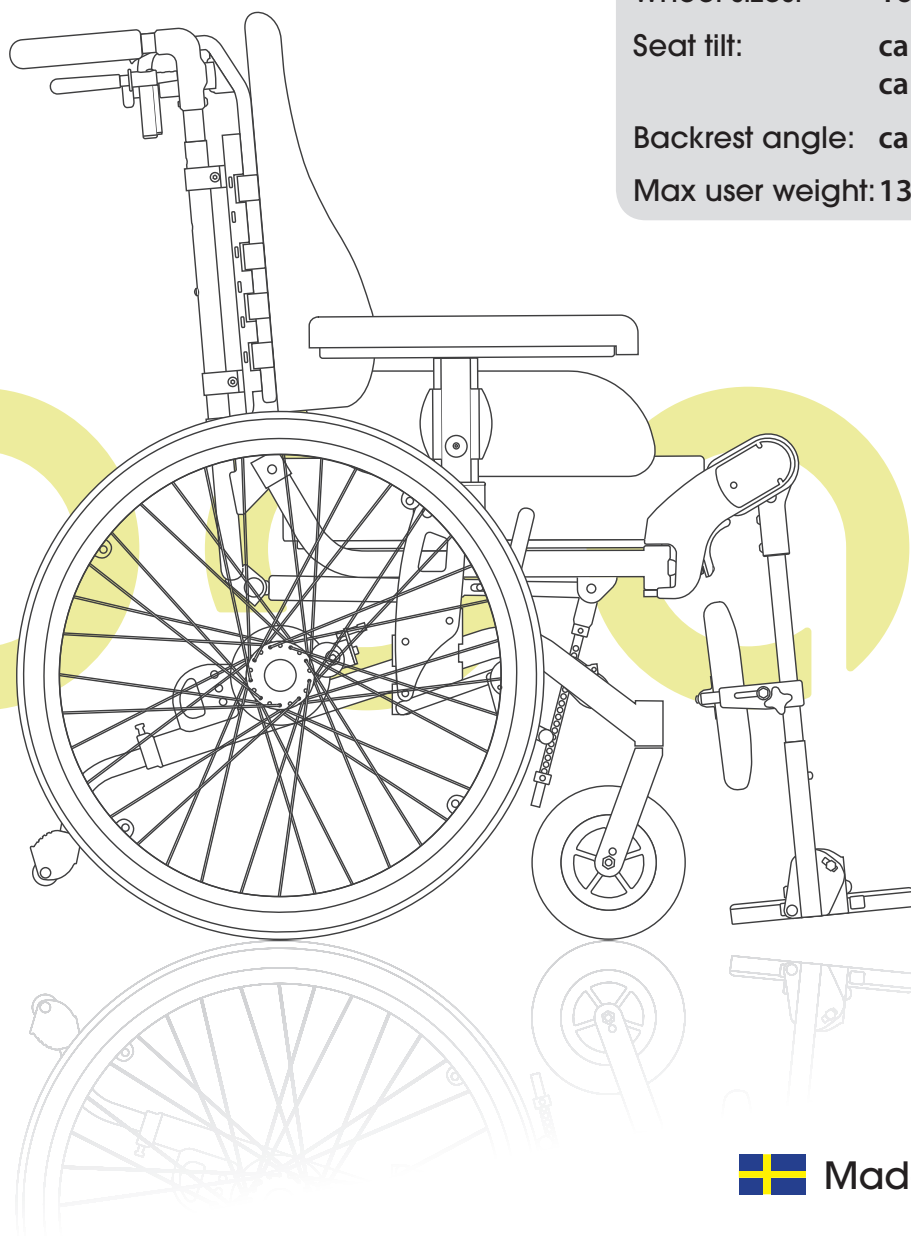

Made in Sweden

## Important features:

### ▶ GOOD HANDLING

Extremely short distance between the caster & wheel gives very good handling for both users & carers. The driving wheels are also fully in line with the castors allowing the chair to run easily on narrow ramps.

### ▶ SEAT TILTING PIVOT IS LOCATED NEAR USER CG

Provides better control of tilting as it requires less force to tilt, especially valuable when tilting heavy users. The tilting pivot position offers the opportunity for expanded tilting in special cases.

### ▶ HIGH PATIENT SECURITY

Scales for different settings

The good seating ergonomics is combined with a number of scales to make it easier for carers and relatives to ensure that the user is sitting in the right position. The communication between prescribers and carer / user / family is greatly simplified by this.

### ▶ CONTROL PLACEMENT

Tilt and back angle controls are positioned so that the carer can maintain his/her grip around the push bar when using these functions. It gives greater control of the movement!

### ▶ FLEXI BACK WITH VELCRO-STRAPS

Flexi back with adjustable Velcro-straps is standard for HD Balance. As an option it is possible to choose increased length and/or width.

### ▶ LEG SUPPORT WITH EXTREMELY GOOD ANGLE ADJUSTMENT

The leg supports can be angled between 90 and 180 degrees by default, without the casters need to be replaced with smaller ones. It is also very easy to remove or swing away the leg support unit. As an option there are leg supports available in different fixed angles.

### ▶ ROBUST AND FLEXIBLE ARMREST

The armrest is mounted in one anchor point which makes it easy to take off and put in place. Setting options in depth and height is generous. The armrest also features a mounting rail where various accessory brackets can be mounted, for example table mounts.

### ▶ A WIDE AND FLEXIBLE RANGE OF ACCESSORIES

A variety of accessories are available such as alternative leg and calf support, trunk support, shoulder support, tables etc..

### ▶ INTEGRATED TRANSPORT MOUNTS

### ▶ THE WHEELCHAIR HAS A CLEAN UNDERBODY AND HIGH GROUND CLEARANCE

